



### PRODUCT INFORMATION

**Product Name:**

PROCEEDx HPV 18Positive Sample

**Catalogue Number:**

VP-62-18

**Product Description:**

This product contains inactivated HPV in storage buffer.

### SPECIFICATIONS

**Strains:**

HPV type 18 wt 4

**Storage:**

Store this preparation at 2 °C to 8 °C.

**Shelf Life:**

24 weeks from manufacturing date.

**Volume:**

1 mL

### USE

**Instructions For Use:**

Use PROCEEDx product as per sample handling instructions specified by device or instrument manufacturer.

**Intended Use:**

The product is to be used **for nucleic acid** test instrument (device) verification/validation, manufacturing QMS, educational and R&D purposes.

**NOT FOR USE IN DIAGNOSTIC PROCEDURES. FOR RESEARCH USE ONLY.**

### ASSITANCE

If you have any questions regarding the production, testing or use of this antigen, please send them together with any relevant data to Microbix Technical Services, at [customer.service@microbix.com](mailto:customer.service@microbix.com) or fax 905-361-8911. Your complete satisfaction with the performance of this product is important to us.

### HAZARDS

The reagent has been inactivated and should contain no infectious material: however, it should be handled as if potentially infectious. Generally accepted good laboratory practices appropriate to biological reagents should be employed. The buffer solution contains methanol (CAS no. 67-56-1) as a preservative. Please refer to the Safety Data Sheet (SDS).

### SPECIAL HANDLING



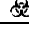
It is recommended that the product is tested according to kit Manufacturer's instructions for testing of clinical samples. Before use, the product should be thoroughly mixed by vortexing for approximately 30 seconds. In protocols for the assessment of samples collected by swab procedures, the amount of liquid should be pipetted into the sample transport tube or a cartridge as per test manufacturers' instructions.

## QUALITY CONTROL INFORMATION

qPCR Test Kit	Target	Result
Cepheid Xpert HPV	HPV-16	Negative
Cepheid Xpert HPV	HPV-18, HPV-45	Positive
Cepheid Xpert HPV	<i>Chlamydia trachomatis</i>	Negative
Cepheid Xpert HPV	Human DNA	Positive
Cepheid Xpert CT/NG	<i>Neisseria gonorrhoeae</i>	Negative

## EXPECTED RESULTS

Testing of the product in PCR-based assays for HPV-18 should result in a positive finding for HPV-18, including such assays where sample adequacy is confirmed for the determination of a positive result.

REF	Catalog Number
LOT	Lot Number
	Temperature Limitation
	Expiration Date
	Biological Risk

